


STATIV[®]

STATIV[®]

*All suture anchor
Optimal bone preservation with
maximum fixation*

STATIV All Suture Anchor

- 1.5mm, 1.8mm, 2.5mm
- Variety of suture/tape combinations
- Minimally invasive
- Ligamentoplasty surgeries

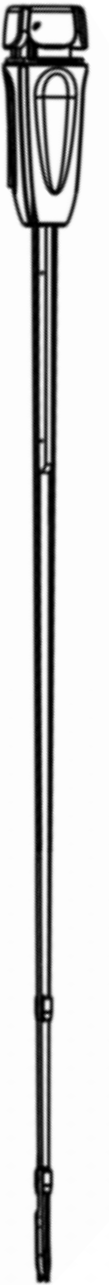


T: 012 807 5172
www.ortho-design.com
info@ortho-design.co.za

ORTHO-DESIGN
ORTHOPAEDIC INNOVATION


Features and Benefits

- A soft anchor with patented design provides benefits of bone preservation while giving fixation benefits of traditional anchors.
- Anchor is loaded on minimalistic designed inserter device with bead to give assured centralized deployment into the pilot hole.
- Small anchor and pilot-hole size allows for multiple anchor placement in difficult insertion areas.
- The deployed unique tri-pod bunched up pattern doubles in size and gives a superior pull out strength.
- Reusable instruments ensure reduced cost and less waste.

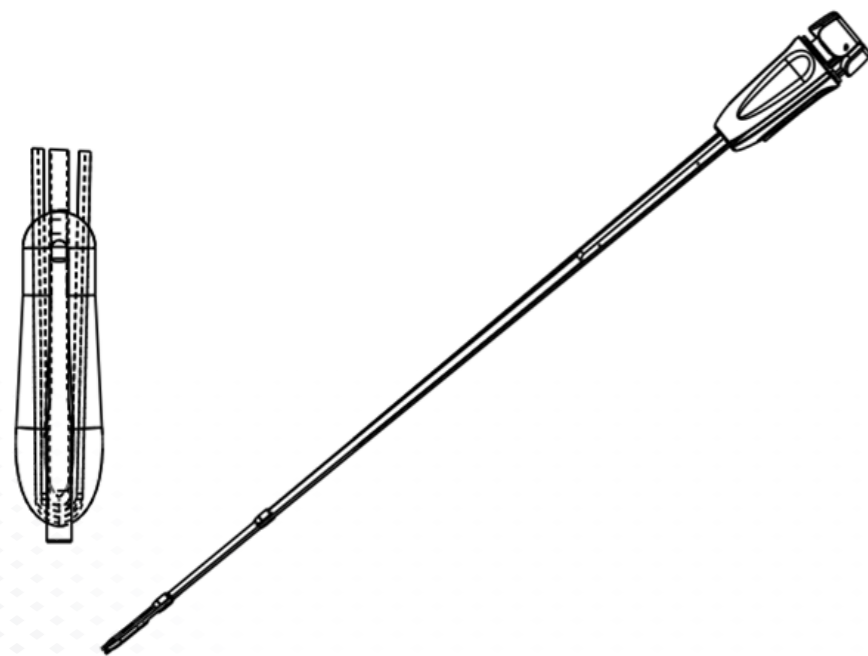

Materials & Specifications

- Anchor Material: Non-absorbable suture anchor-braided, UHMWPE
- Suture Material: Braided, UHMWPE non-absorbable #2 suture
- Tape Material: Braided, UHMWPE non-absorbable 1.5mm flat tape
- Handle: ABS plastic
- Introducer: Stainless Steel

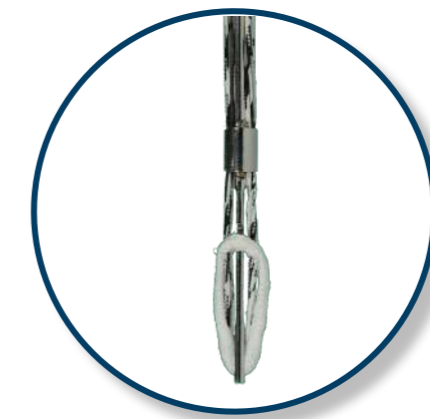

Extensive Application Range

- The STATIV All Suture Anchors are intended to be used for suture or tissue fixation in the foot, ankle, knee, hand, wrist, elbow, shoulder and hip.
- Highly diverse range provides extensive product indications.

Available Options:

Size	1.5mm	1.8mm	2.5mm
Suture Loading	Single Suture	Double Suture	Double/Triple Suture
Tape Loading	Single Tape	Double Tape	Double Tape