


VersaQuad™

Implant System
2 x VersaTap™ Suture Anchors
2 x VersaLat™ Knotless Anchors

- 2 x VersaTap™ Anchors**
- two #2 UHMWPE sutures
 - single-use inserter

- 2 x 6mm VersaLat™ Anchors**
- one #2 UHMWPE suture
 - single-use inserter

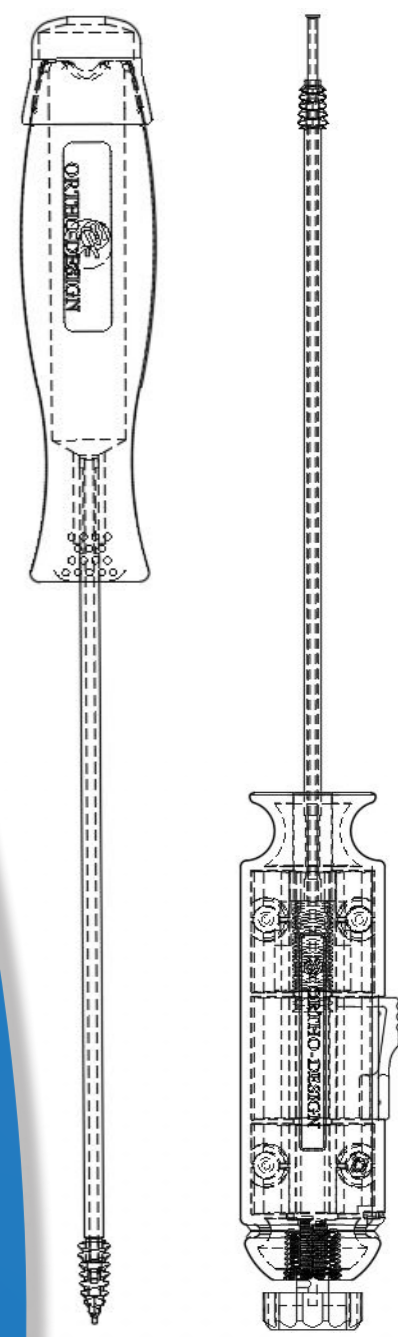

T: 012 807 5172
www.ortho-design.com
info@ortho-design.co.za



Features and Benefits

VersaTap™ Features and Benefits:

- Self-tapping titanium anchor tip integrated with a unique cutting edge that eliminates the need for a pre-tapped hole.
- A multi-start thread leads to effortless rigid insertion.
- Tapering thread design for improved fixation between the anchor and the cortical bone.

VersaLat™ Features and Benefits:

- A one-handed technique promotes ease-of-use and decreases operation time.
- Allows knotless fixation of a large number of sutures due to the characteristic suture loop.
- Unique locking mechanism for effective & efficient suture tensioning.
- Eliminates knots associated with tissue irritation.
- A retention suture that allows for the fixation of any remaining ligament fragments.



Materials & Specifications

VersaTap™ Materials & Specifications:

- Anchor size: 6.75mm x 17.5mm
- Sutures: Two #2 Force Fibre UHMWPE braided suture (non-absorbable)

VersaLat™ Materials & Specifications:

- Anchor Size: 6 x 12mm
- Retention Suture Material: One #2 Force Fibre UHMWPE braided suture (non-absorbable)

Extensive Application Range

- The VersaQuad™ implant system offers a more affordable double-row rotator cuff repair solution for patients. Note: one of the VersaLat™ anchors can be used to perform a biceps tenodesis.
- The VersaTap™ is used as a medial anchor row in double-row rotator cuff repair surgery.
- The 6mm VersaLat™ to be used as a knotless anchor in single and double-row rotator cuff repair surgery.

