




VersaTap™

The future of effortless, rigid fixation

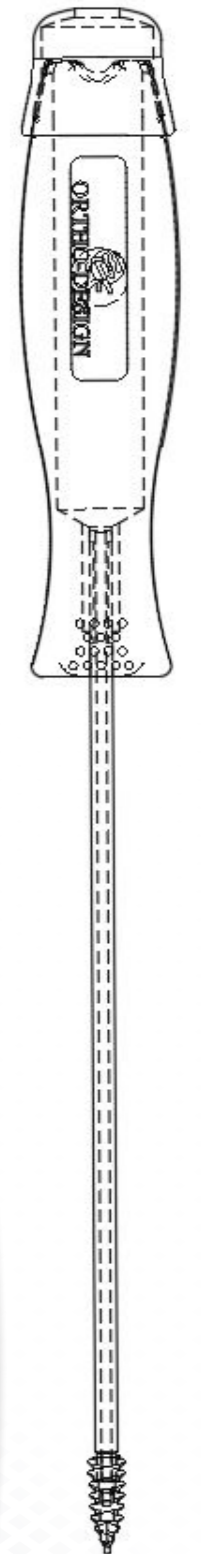

Multi-compound VersaTap™

- two #2 UHMWPE sutures
- single-use inserter

Multi-compound VersaTap™

- three #2 UHMWPE sutures
- single-use inserter

Multi-compound VersaTap™

- two #2 UHMWPE sutures
- single-use inserter
- half round taper point needles

Multi-compound VersaTap™ with Tapes

- two #2 UHMWPE tapes
- one #2 UHMWPE suture + one #2 UHMPWE tape
- single-use inserter




T: 012 807 5172
www.ortho-design.com
info@ortho-design.co.za

ORTHO-DESIGN
ORTHOPAEDIC INNOVATION

Features and Benefits

- Self-tapping titanium anchor tip integrated with a unique cutting edge that eliminates the need for a pre-tapped hole
- A multi-start thread that leads to effortless rigid insertion
- Tapering thread design for improved fixation between the anchor and the cortical bone
- A PEEK anchor body that results in minimal post-operative artifacts
- PEEK Body allows for clear imaging of the ligament bone interface
- Reduced environmental impact, wastage and cost when using our reusable inserter, sustainable single-use handle material

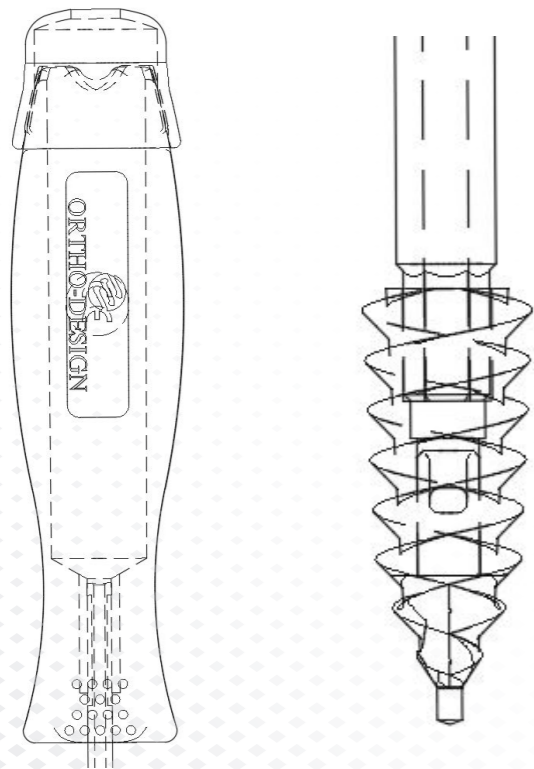

Materials & Specifications

- Anchor Body Material: 100% Medical grade Polyether Ether Ketone (PEEK)
- Anchor Tip Material: Medical grade titanium tip (Ti-6Al-4V)
- Anchor size: 6.75mm x 17.5mm
- Sutures: Two #2 Force Fibre UHMWPE braided suture (non-absorbable)



Extensive Application Range

- Used in single-row rotator cuff repair surgery
- Used as a medial anchor row in double-row rotator cuff repair surgery
- Various anchor sizes and suture options to accommodate patient-specific needs
- Introducer options available to accommodate open and arthroscopic repair
- Well suited for all general surgical solutions

