


VersaTi™ BROCHURE

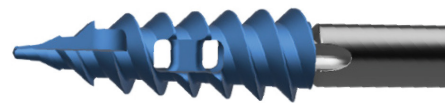

VERSATI™

The VersaTi™ follows the same design methodology as the VersaTAP™. This anchor range is available in certified medical grade titanium and can be supplied with a drill in applications where thicker cortical bone is present.

FEATURES AND BENEFITS

The self-tapping titanium tip with a unique cutting edge eliminates the need for a pre-tapped hole

A multi-start thread leads to effortless insertion

A tapering thread design improves fixation between the anchor and the bone

A unique retractable silicone handle for effortless and efficient open repair surgery



Available Sizes:

Item Code	Description	Sizes
ATP001	VersaTi, Titanium Suture Anchor	5.5mm x 17mm, W/Needles
ATP001-CH	VersaTi, Titanium Suture Anchor	5.5 x 17mm
ATP004	VersaTi, Titanium Suture Anchor	3 x 11mm, W/Needles

References

ATP001	VersaTi, Titanium Suture Anchor, 5.5 x 17mm, w/Needles
ATP001-CH	VersaTi, Titanium Suture Anchor, 5.5 x 17mm
ATP004	VersaTi, Titanium Suture Anchor, 3 x 11mm, w/Needles

MATERIALS & SPECIFICATIONS

Anchor Material: Certified ISO 5832-3 Medical Grade Titanium (Ti-6Al-4V)

Anchor size: Various sizes available

Suture Material: Two #2 Force Fibre UHMWPE braided sutures (non-absorbable)

Needles: Four M06, half-round taper point needles (alternative options available)

Introducer: Single-use, open-handle introducer composing of medical grade stainless steel and silicone

EXTENSIVE APPLICATION RANGE

Used in single-row rotator cuff repair surgery

Used as a medial row anchor in double-row rotator cuff repair

The anchor range can be used in shoulder, elbow, wrist, knee and ankle surgery

Anchor size and shape variations to ensure optimal fixation for large and small-joint repair

Well suited for all general surgical solutions