

VersaTap™

SUTURE ANCHOR



VERSATAP™

The VersaTAP™ is a self-tapping PEEK anchor with a titanium tip. This anchor is ingeniously designed to combine the advantages of both PEEK and Titanium anchors. The titanium tip allows for self-tapping into the bone, whilst the PEEK body allows for minimal post-operating artifacts.

FEATURES AND BENEFITS

Self-tapping titanium anchor tip integrated with a unique cutting edge that eliminates the need for a pre-tapped hole


A multi-start thread that leads to effortless rigid insertion

Tapering thread design for improved fixation between the anchor and the cortical bone

A PEEK anchor body that results in minimal post-operative artifacts

PEEK Body allows for clear imaging of the ligament bone interface

Reduced environmental impact, wastage and cost when using our reusable inserter, sustainable single-use handle material

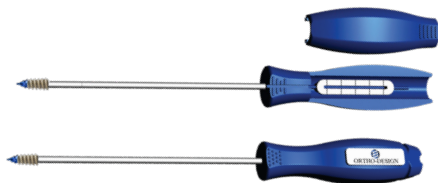

MATERIALS & SPECIFICATIONS

- Anchor Body Material: 100% Medical grade Polyether Ether Ketone (PEEK)
- Anchor Tip Material: Medical grade titanium tip (Ti-6Al-4V)
- Anchor size: 6mm x 17mm
- Sutures: Two #2 Force Fibre UHMWPE braided suture (non-absorbable)

EXTENSIVE APPLICATION RANGE

- Used in single-row rotator cuff repair surgery
- Used as a medial anchor row in double-row rotator cuff repair surgery
- Various anchor sizes and suture options to accommodate patient-specific needs
- Introducer options available to accommodate open and arthroscopic repair
- Well suited for all general surgical solutions

Available Configurations

ADP001	VersaTap™, PEEK/Ti Suture Anchor	6 x 17mm	w/2 Sutures
ADP001-T1	VersaTap™, PEEK/Ti Suture Anchor	6 x 17mm	w/1 Tape, 1 Suture
ADP001-T2	VersaTap™, PEEK/Ti Suture Anchor	6 x 17mm	w/2 Tapes
ADP002	VersaTap™, PEEK/Ti Suture Anchor	6 x 17mm	w/3 Sutures
ADP005	VersaTap™, PEEK/Ti Suture Anchor	6 x 17mm	w/2 Sutures, 4 Needles
ADP001-LT1	VersaTap™, PEEK/Ti Suture Anchor	6 x 17mm	w/1 Tape
ADP055	VersaTap™, PEEK/Ti Suture Anchor	5.5 x 17.5mm	w/2 Sutures

ORTHO-DESIGN

ORTH-03-08-R-VT-23 (REV 2)