


VersaLat™

KNOTLESS ANCHOR

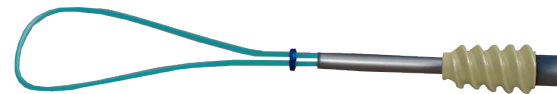

VERSALAT

The VersaLAT™ is an innovative lateral knotless anchor. This knotless screw-in anchor can be combined with any medial row suture anchor or as a stand-alone lateral anchor.

FEATURES AND BENEFITS

A one-handed technique promotes ease-of-use and decreases operation time.

Allows knotless fixation of a large number of sutures due to the adjustable suture loop.

Unique locking mechanism for effective & efficient suture tensioning.

Eliminates knots associated with tissue irritation.

A retention suture that allows for the fixation of any remaining ligament fragments.

Engineered counter thread allows for independent rotation of the inner shaft upon anchor insertion.

Available Sizes:

Item Code	Description	Sizes
ADP021-P	VersaLat, PEEK Knotless Anchor	4.75 x 13mm
ADP021-TI	VersaLat, Titanium Knotless Anchor	4.75 x 13mm
ADP008	VersaLat, PEEK Knotless Anchor	5.5 x 12mm
ADP016	VersaLat, PEEK Knotless Anchor	6.5 x 12mm
ADP019	VersaLat, PEEK Knotless Anchor	7.5 x 13mm

MATERIALS & SPECIFICATIONS

Anchor Material: 100% Medical grade Polyether Ether Ketone (PEEK)

Retention Suture Material: One #2 Force Fibre UHMWPE braided suture (non-absorbable)

Introducer: Single-use, cannulated introducer composing of medical grade stainless steel and silicone

EXTENSIVE APPLICATION RANGE

The 4.75mm and 6mm VersaLat™ to be used as a knotless anchor in single and double-row rotator cuff repair surgery.

The 7mm and 8mm VersaLat™ to be used for biceps tenodesis repair.

Can be used as a general interference screw or internal brace.

References

ADP021-P	VersaLat, PEEK Knotless Anchor, 4.75 x 13mm
ADP021-TI	VersaLat, Titanium Knotless Anchor, 4.75 x 13mm
ADP008	VersaLat, PEEK Knotless Anchor, 5.5 x 12mm
ADP016	VersaLat, PEEK Knotless Anchor, 6.5 x 12mm
ADP019	VersaLat, PEEK Knotless Anchor, 7.5 x 13mm