


DueLock™ Knotted Push-In Suture Anchor

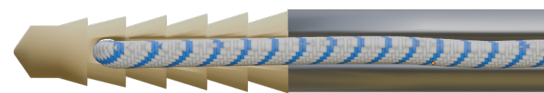

DueLock™

PEEK Push-in Suture Anchor

(Class IIb)

The 2.9mm PEEK DueLock suture anchor has a tapered thread profile uniquely designed to increase pull-out strength and minimise the incision footprint. The associated instrumentation allows for accurate and effortless incision onto the Glenoid rim.

The anchor comes preloaded with a multicoloured ForceFibre suture for improved suture management.

ADVANTAGES:



Knotted Slot:

- Provides the surgeon with a knotted fixation method
- Increased pull-out strength and a reduced footprint
- A tapered thread design to maximize cortical compression

Introducer:

- A range of introducer lengths are available to accommodate small and large-joint repair

PRODUCT ANCILLARY:

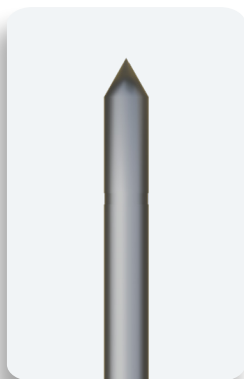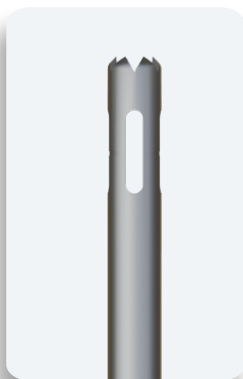

Available sizes:

2.9mm x 10mm

Varying Insert Lengths

INDICATIONS:

The DueLock™ is intended for Bankart, Labral and SLAP repair. The DueLock™ can also be used in hip labral repair and general soft tissue fixation applications.

COMPOSITION:

Anchor: PEEK

Insert: Stainless Steel

Suture: UHMWPE braided, non-absorbable

