


# EcoLat™ SUTURE ANCHOR

---

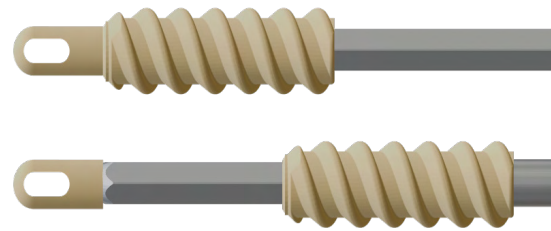


## EcoLat™ Suture Anchor

The EcoLat™ Knotless Suture Anchor is a re-usable system specifically developed for knotless fixation of soft tissue and/or external sutures – applicable to a number of ligamentoplasty indications. This novel implant design, which makes use of a reusable introducer, was designed from the ground up with a dual focus on quality and reduced cost.

### FEATURES AND BENEFITS

- Leading suture eyelet for fine-tuned tensioning
- Can handle up to 8 suture/tape tails
- Square driver pocket maximizes insertion strength
- Reduce implant cost by up to 50-75%
- Reusable instrumentation designed for use with both anchors sizes reduces cost and waste
- In-surgery assembly allows for modularity of suture/tape configurations and reduces inventory burden.
- Vented design allows for blood flow to repair site to aid healing response.

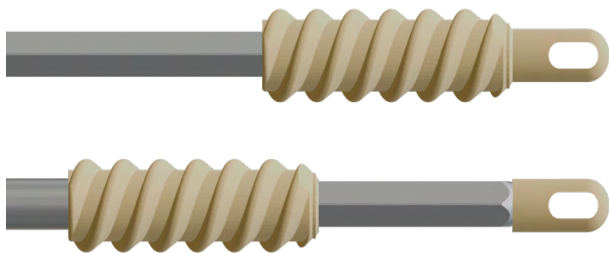

#### Available Configurations:

ECO-LAT-475	EcoLat™, PEEK Knotless Suture Anchor, 4.75 x 22mm
ECO-LAT-55	EcoLat™, PEEK Knotless Suture Anchor, 5.5 x 22mm

### INDICATIONS:

- Rotator Cuff Repair
- Biceps Repairs and Tenodesis
- Ligament and Tendon Repairs

### COMPOSITION:

4.75mm and 5.5mm sizes

Pre-loaded with 1 or 2 HS Fiber UHMWPE USP2

Loaded with 1 or 2 HS Fiber UHMWPE TAPE 1.6mm

#### Anchor Body:

100% PEEK (PolyEther Ether Kentone)

#### Suture Material:

UHMWPE braided, non-absorbable

Solvay's Zeniva PEEK material provides excellent biocompatibility and structural integrity

### SUTURE VARIATIONS:

- HS Fiber UHMWPE White/Blue USP2
- HS Fiber UHMWPE White/Blue TAPE 1.6mm